
Quick guide

Preparations

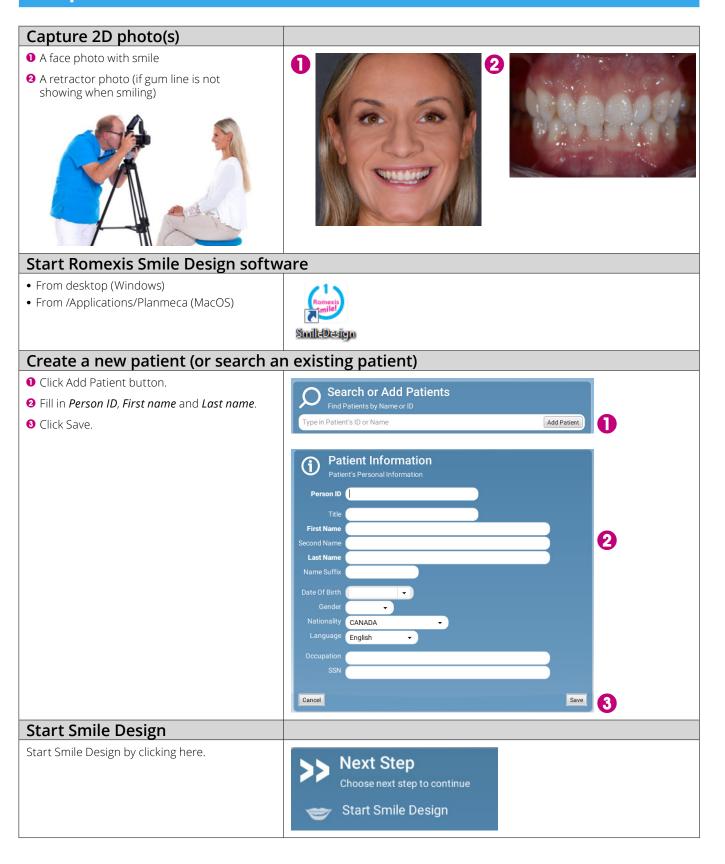

Smile designing

Import face photo

Click Import button and browse the image file.





Align and crop photo


- 1 Click the Align and crop button.
- 2 Drag the blue dots onto two reference points, e.g. the pupils.
- Double-click on the green control point in the middle to align.
- 4 For cropping drag the borders of the image.
- **6** Confirm the crop by double right-clicking.





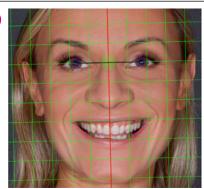

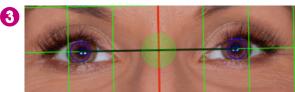

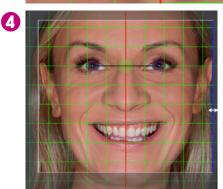

Calibrate image

- Measure the length of e.g the patient's incisor.
- 2 Click the Calibrate button.
- 3 Draw a calibration line onto the image.
- 4 Type the width of the incisor.

















Add teeth design

- Add a teeth design to the image by clicking the Teeth button.
- Open a suitable shape from the teeth library by clicking the black triangle on the Teeth button.











Scale, translate, bend and rotate teeth as a group

- To scale, translate, bend and rotate the teeth as a whole use the blue and red control points.
- 2 To resize all teeth simultaneously drag the blue rectangle around the teeth.
- **3** To delete the teeth design click the Delete button in Smile Design toolbar.








Rotate, scale and reshape individual tooth


- Activate the tooth by clicking on it.
- Rotate the tooth from the red control point.
- Resize the tooth by dragging the blue rectangle.
- Reshape the tooth by dragging from the blue control points.
- To add a control point double-click on the tooth outline, to remove a control point double-click on it.

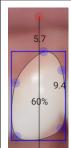

Mirror, hide and reset tooth

Using the rightclick menu you can:

- Mirror one tooth to the opposite side
- Mirror all teeth to the opposite side
- Decide how may teeth to include in the design – Hide/show tooth/all teeth
- Show Golden proportions ratios
- Save teeth designs to library
- Select a colour for each tooth separately

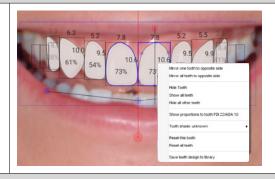

Define smile area

- 1 Click the Smile Area button.
- 2 Drag the line to translate or the control points to reshape.

Double-click to add/remove points.







Select teeth shade


- Select shade from standard VITA classical shades and bleached shades.
- Select shade from the image using the Colour picker tool.
- 3 Adjust shade tone.





Visualise teeth

- Adjust transparency.
- 2 Adjust photorealistic texture.
- 1 Show/hide teeth design.
- 2 Show/hide smile area.

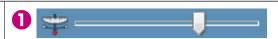



View before and after images in full screen view






Using advanced tools

Clone brush

- Click the Clone brush button.
- 2 Select the clone mode and the suitable settings.
- Olone the image area by right-clicking on a source area and clone it on desired area by dragging with the mouse.
- Select Warp Brush to stretch tissue by dragging with the mouse.
- **3** Select Mirror to clone in mirroring mode.

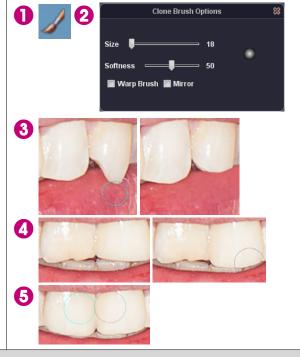


Align retractor image

- 1 Import intraoral image by clicking this button.
- Select retractor image from folder. The images open up in matching dialogue.
- Opening the state of the sta
- Check the positioning by adjusting transparency. Click OK.
- 6 Adjust teeth and intraoral image




Add drawings and measurements

- Add a generic implant, text, lines and measurements to the image using these tools.
- 2 To select annotation click the selection button and then click the annotation. Selected annotations can be modified and deleted.

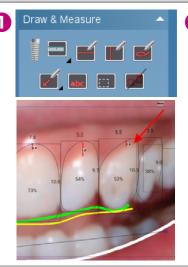

Export options

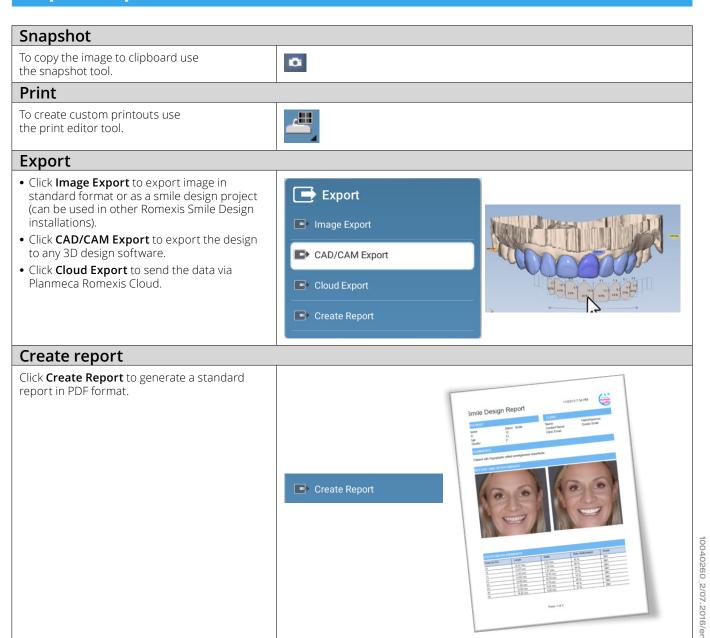

To read the full Planmeca Romexis SmileDesign manual visit:

http://manualkit.planmeca.com/#

To watch tutorial videos:

Click this button at the top right corner of the screen or browse to

www.youtube.com > search for "Planmeca Smile Design" playlist https://www.youtube.com/watch?v=7N1ub5oe98w&list=PLPL3t-3CNS9U7Oo2JRskt4t5GFUvY-yax